Notice

This product is intended for household use and should not be used for industrial or medical purposes. During use, you should make sure to protect your health, consult your doctor before using this product.

- Please follow the instructions to operate the product, don't use other parts without the Company's recognition.
- Do not use the machine if you having pregnancy, a wound, skin inflammation, skin diseases, heart diseases, etc.
- Don't use it within one hour after a meal.
- Children or disable people should have adult supervision when using it.
- Do not expose the machine at high temperatures or wet environment.
- Do not use the machine in damp places, such as the bathroom.
- Please cut off the power during storage.
- Do not pull or pull the power line to the machine.
- Do not let water, oil or liquid seep into the machine. (Massage head can be cleaned by water)
- Do not use corrosive the plastics and metal detergent to clean the machine.

Use the cloth cover, protect the hair.

Unuse the cloth cover, the hair be dangerous.

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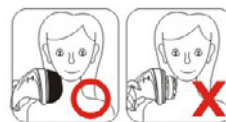
Main Function

- Light weight body design make easy to use it.
- Quiet design make comfortable.
- With stepless variable speed design, the strength of massage is adjustable.
- Powerful and truly deep massage with high-torque motor.
- Add three massage heads (smooth, wave, ball), both decompression, massage, essential oils pushed into grease and other functions. And add protective cloth cover so that prevent hair involved in the machine when you massage the neck.
- Demolition massage head design that allows you to assembly massage head easily, and more convenient to clean massage head.

Specification

Product size: approximately 199 x 119 x 151 mm

Main machine material: ABS/TRP plastic

Standard Accessories: motor, smooth type massage head, wave type massage head, ball type massage, and protective cloth cover.

Multiple sets of massage heads, a set with multi-usage

- 1. Smooth massage head:** Recommend to massage the neck to relieve stress.
- 2. Wave type massage head:** acupuncture points and deep tissue massage, strengthening to improve the blood circulation and body functions, like physical therapy professional masseurs to enjoy.
- 3. Ball type massage head:** oil massage, powerful import, in depth activation of absorption efficiency, **do not install the protective cloth cover.**
- 4. Protective cloth cover:** prevent hair involved in the machine from the head, **not suit to ball type head massage.**